
WESTMED

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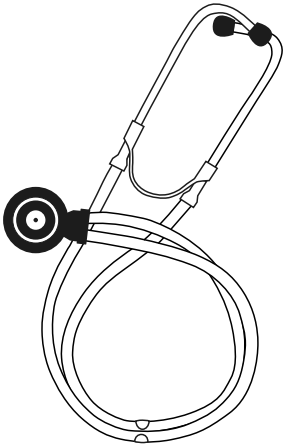EN

by manufacturer

Please read these instructions for use carefully and completely before use and observe the care instructions.

This medical device is exclusively manufactured in accordance with the regulations in force, relating to medical devices and is subject to quality controls.

Instructions for use

Please read carefully the user manual before using the product and keep it safe for further reference.

Intended medical use

The Rappaport-Stethoscope is a professional equipment for aid in the diagnosis. It is an acoustic instrument used for auscultation to listen to sounds within the human body in order to make a medical diagnosis: for example during heart, lungs and other organic noises auscultation as well as measuring the blood pressure.

Indications


The Rappaport-Stethoscope is a medical device Class I, for aid in the diagnosis, intended for health care professionals, trained for its use. It is used to listen to sounds within the body for gathering information about patient condition and to measure blood pressure for all types of patients: adults, children and newborns.

Contraindications

- Do not used on injured skin.
- Do not used in case of instruments material allergies.

Instructions for Use

Fitting and adjusting the binaural tubes.

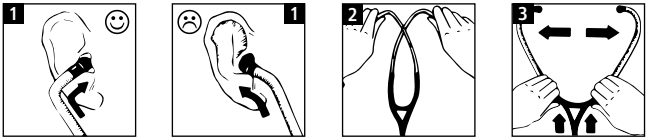The binaural tubes are designed to automatically adjust to your ears listening angle. Ensure that the branches of the binaural tubes are pointing forwards, thus matching the direction of the auditory canal [fig. 1].
The binaural tube spring can be adjusted for improved comfort and better adjustment to your morphology.
To increase the spring tension, grasp the branches of the binaural tubes at the level of the earpieces and cross them over [fig. 2].
To reduce the spring tension, grasp the branches at the base of the binaurals and gradually move them apart to complete extension [fig. 3].

Selecting chestpiece

To select the desired chestpiece, gently tap this latter‘s membrane with the tip of your finger. If there is no sound, rotate the chestpiece through 180 ° until a click is heard.

Using the big chestpiece

The big chestpiece allows high frequency sounds to be heard.

Using the small bell chestpiece

The small bell chestpiece allows listening to low frequencies sounds.

Material	Anodized Aluminium
Weight	37,5 g
Diameter Diaphragm	47 mm
Diameter Bell	30 mm
Diaphragm	
Material	PVC
Diameter	47 mm
Non-chilling ring	
Material	Thermo-propylen rubber
Diameter Diaphragm	47,3 mm
Diameter Bell	37,5 mm
Binaurals and other parts	
Tubing length	60 cm
Tubing material	PVC
Binaurals material	Copper
Spring material	Stainless steel
Ear tips material	Metal
Total length	76-77 cm
Accessories included	1 pair of ear-tips

Warnings

- Do not immerse the stethoscope in water/liquid.
- Do not use if the packaging is damaged.
- Do not sterilize.

Side effects

None.

Materiovicilance

All serious incidents related to the use of a medical device should be notified to the manufacturer and the responsible authority of the Member State where the user and/or the patient live.

Storage and handling conditions

- Keep the stethoscope away from excessive heat or cold.
- Avoid exposing stethoscope to direct sunlight.

Handling conditions	Temperature	+10 °C to +40 °C
	Humidity	20% to 85%
Storage conditions	Temperature	-20 °C to +70 °C
	Humidity	20% to 85%

Upkeeping and cleaning recommendations

Regular cleaning of your stethoscope and in particular of the earpieces. ensures optimum acoustic performance.

- Clean the stethoscope with a cloth soaked in 70% isopropyl alcohol or moderately soapy water.
- Unscrew the earpieces before cleaning them thoroughly by soaking, then drying them Ensure that no dirt remains inside.

Disposal

Dispose of the device at the specific collection point according to local regulation.

Guarantee

The Rappaport-Stethoscope is guaranteed against manufacturing defects (parts and labour) under normal conditions of use, from the date of purchase (invoice date).

- Stethoscope: 1 year.

The guarantee does not apply in the event of use non-compliant with instructions. improper use. modifications made by yourself or by an unauthorised third party. accidental damage and normal wear. During the guarantee period, any stethoscope parts that present manufacturing defects shall be replaced and/or repaired free of charge. after returning the device to Spengler's after-sales service. In the event of a guarantee claim or of repairs, please return the complete stethoscope in its original packaging, accompanied by the dated invoice.

This information is provided without guarantee and is subject to change.

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Fleischhacker GmbH & Co. KG
Ltd No. 8 Tongxing Road Economic & Technological
Shanghai International Holding Corp. GmbH (Europe)
MedNet SWISS GmbH